



**RELAX AND RELAX 2. RELAXATION CHAIRS FOR RETIREMENT HOMES,
NURSING HOMES, HOSPITALS AND CLINICS.**

 **GREINER**

ambient
care.



GREINER CHAIRS, RECLINERS AND RELAXATION CHAIRS FOR GERIATRIC CARE AND MOBILISATION.

A good seat makes a lot of things easier – pure relaxation is quality of life. Sitting and lying comfortably for shorter and longer periods during an active daily routine contributes greatly to well-being and satisfaction for older people in care. At the same time, satisfied residents and patients enormously reduce the load on nursing staff. People need care for various reasons: after a long time in hospital; when they find it difficult to cope alone; when their mental powers diminish. In many areas of modern care facilities, GREINER recliners and relaxation chairs create a new form of home and nursing culture in an appealing ambience. The spectrum ranges from day care via assisted living to geriatric care – from rest and relaxation to mobilisation and transport.



CareLine.
Care and mobilisation chair.
See website.

RELAX AND RELAX 2.

Design, comfort and safety are what give nursing chairs and relaxation chairs their feel-good factor. With comfortable upholstery, an appealing choice of colour for the covers, easy adjustment of the sitting and lying position as well as accessories such as an attachable wooden table, the beautifully designed GREINER seating furniture is a defining element in rooms and communal areas. RELAX and RELAX 2 were developed specifically for the residential and care needs of elderly care and nursing institutions but are also frequently used in hospitals. Comfort for the residents while seated or lying down is optimally combined with requirements on the carer's side for stability, safety and hygiene.

Recliners and relaxation chairs RELAX and RELAX 2

- ◆ A wide choice due to models with manual or electric adjustment of the reclining position, stationary and mobile versions as a variety of accessories
- ◆ A high degree of comfort thanks to smoothly adjustable backrests and legrests, right up to the resting position, an ergonomic shape, cushioning using viscous foam and optional premium cover fabrics in many colours
- ◆ A pleasant living atmosphere thanks to the combination of natural beech wood with attractively designed and hygienic easy-care surfaces – coated with polyurethane-based varnish
- ◆ Sturdy armrests with a curved shape ensure safety and make it easier to get up from the resting position
- ◆ Easy to clean due to removable seat cushions and covers



RELAX. RELAXATION CHAIR FOR RETIREMENT HOMES, NURSING HOMES, HOSPITALS AND CLINICS.

The RELAX relaxation chair ensures all-round comfort. The backrest and legrest are continuously and synchronously adjustable to about 50° from the neutral position, with the option of electrical or manual adjustment via a gas spring. In addition, it is possible to separately adjust the leg section by up to 90°. In the relaxed position, the seat cushion rises at the front by 10°, thus preventing the seated person from slipping. The rocking function is particularly relaxing and is easily operated using the adjustment lever under the armrest. The seat cushion and cover are removable. The shape of the armrest provides safety and comfort in reclining position. RELAX is available as a fixed version and as a mobile version on double castors. These can be locked on the rear side.



Mobility with RELAX in nursing



RELAX with attached patient table



RELAX without side panels in resting position



A feel-good atmosphere with RELAX in nursing homes or in hospitals



Head cushion with side bolsters (accessory)



Versatile colour variants in the cover fabrics compact and premium

relax – Accessories

• Cover materials	Premium material according to our colour chart, also for head cushion and armrests
• Head cushion	With/without side bolsters
• Armrests	Padded (pair)
• Side elements	Closed, for lateral support (pair)
• Support cushion	Secures the back (pair)
• Push bar	Beech
• Patient table	Beech plywood, detachable

Dimensions and technical data

• Seat width	48 cm
• Seat depth	49 cm
• Total depth in basic position	81 cm
• Total depth in lying position	160 cm
• Total height in basic position	119 cm
• Total height in lying position	90 cm
• Lying surface	170 cm
• Mounting height	46 cm
• Maximum patient weight	125 kg
• Weight	30 – 32 kg
• Total width	63 cm
• Roller diameter	75 mm

relax 2.



RELAX 2. PREMIUM RECLINER FOR HOMES AND NURSING.

In the form of an armchair, RELAX 2 impresses with a completely new technology for backrest and legrest adjustment using two gas springs, which allow both synchronous adjustment of the lying surface and separate adjustment of the legrest in every position. The electrically adjustable version is equipped with 2 motors. A powerful battery guarantees the required mobility. Extension of the adjustable backrest into the horizontal lying position, the visco-elastic foam embedded in the upholstery and the ergonomic design offer maximum comfort while seated or lying down. The side panels made of beech wood, polished with polyurethane-based varnish and suitable for wipe disinfection, are attractively stained in the shade Havana. The broad sweep of the armrests beyond the backrest offers complete support of the arms even when lying down and makes it easier to get up from the resting position. 75 mm swivel castors, lockable at the rear, and an ergonomically fitted push handle (accessory) facilitate mobility for the nursing team.



Mobility with RELAX 2 in nursing



RELAX 2 with attached patient table




RELAX 2 without side panels in the resting position



Feel-good atmosphere with RELAX 2 in nursing homes or in hospitals



Ergonomically fitted push bar (accessory)

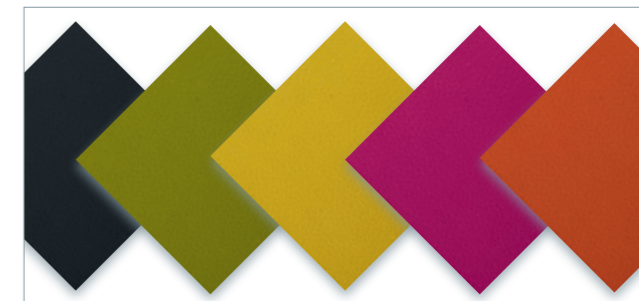

Versatile colour variants in the cover fabrics compact and premium

relax 2 – Accessories

• Cover fabrics	Premium material according to our colour chart, also for head cushion and armrests
• Head cushion	Moveable and detachable
• Armrests	Padded (pair)
• Side panels	Closed, for lateral support (pair)
• Support cushion	Secures the back (pair)
• Push bar	Beech stained Havana
• Patient table	Beech plywood, detachable

Dimensions and technical data

• Seat width	51.5 cm
• Seat depth	48.5 cm
• Total depth in basic position	82 cm
• Total depth in lying position	169 cm
• Total height in basic position	130 cm
• Total height in lying position	87 cm
• Lying surface	180 cm
• Access height	47 cm
• Maximum patient weight	150 kg
• Weight	55 – 65 kg
• Total width	72 cm
• Castor diameter	75 mm



we design quality.

Agreeably health-promoting seating has motivated us since 1922. For decades we have been making high-quality, technically sophisticated couches within the medbest range. Our innovative devices for examination, dialysis, oncology, blood sampling, care and transport are all made in Germany. Our quality commitment is verified by an independently audited quality assurance system: GREINER is certified under DIN EN ISO 9001 and 13485.



Test mark for
quality management system



GREINER

GREINER GmbH
Wettestraße 1 · 74385 Pleidelsheim · Germany
T +49 7144 81120 · F +49 7144 8112 99
info@greiner-gmbh.de · www.greiner-gmbh.de