

sita.



SITA RECLINER CHAIR FOR GENERAL PRACTICE, OFFICE AND CLINICS.

 **GREINER**

sita.



SITA RECLINER CHAIR FOR GENERAL PRACTICE, OFFICE AND CLINICS.

SITA is the chameleon among chairs for clinics and surgeries. Its adaptability to a wide range of applications, supported by versatile accessories, accentuates perfectly its attractive price-performance ratio. Whether fixed or mobile on castors, it is used primarily for blood collection and infusion therapy, for treatment or as a rest and acupuncture chair. The generously dimensioned upholstered surfaces, optionally available in many colours in the super-soft GREINER premium material, the upholstery and the softly-constructed arm rests ensure a high level of patient comfort.


Operation is as simple as it is convenient: when the right-hand operating lever is pushed back, the mechanism unlocks the backrest and can be steplessly moved horizontally. When discharged, it self-aligns thanks to a pneumatic spring. The forward lever position regulates the leg rest, also supported by a pneumatic spring. Perfect ergonomics for patients. And helpful relief for the nursing staff, whose tasks are best supported by a stable base with capacity for 220 kg patient weight, individually steerable 10 cm castors and a handle on the mobile model, as well as accessories such as phlebotomy armrest, infusion rod, patient table and the standard rails for attaching individual utensils.



SITA – perfect for patients and caregivers




Robust base with smooth-running castors

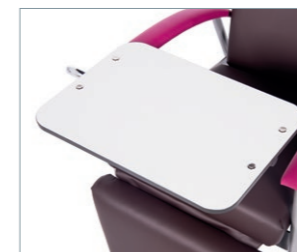

Patient table detachable and movable, at right hand side



Blood collection armrest with ball joint

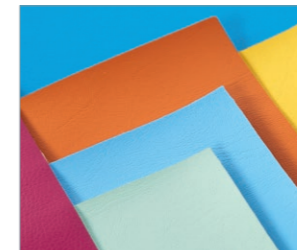

Versatile colour options in the compact and premium upholstery fabrics



SITA in premium material mocha/fuchsia with blood collection armrest on the standard rail and infusion rod in holder



Access position



Seating position



Relax position



Bed position/Flat position

SITA – accessories

- Premium material according to colour chart (upholstery, headrest and padded armrests)
- Headrests in premium material
- Armrests upholstered in premium material
- Phlebotomy armrest in premium material
- Standard rails (pair)
- Phlebotomy armrest and holder for attachment to standard rail with double ball joint, swivable (1 piece)
- Infusion rod with holder
- Patient table detachable and movable, mounted on the right hand side
- Urine bag holder left and right (pair)
- Transparent protective film for leg rest

Dimensions and technical data

• Seat width	50 cm
• Seat depth	49 cm
• Total depth in a sitting position	82 cm
• Total depth in reclined position	185 cm
• Total height (without accessories)	126 cm
• Bed surface	185 x 50 cm
• Seat height	55 cm
• Maximum patient weight	220 kg
• Weight	41 kg
• Total width	63 cm
• Roller diameter	100 mm



we design quality.

Agreeably health-promoting seating has motivated us since 1922. For decades we have been making high-quality, technically sophisticated couches within the medbest range. Our innovative devices for examination, dialysis, oncology, blood sampling, care and transport are all made in Germany. Our quality commitment is verified by an independently audited quality assurance system: GREINER is certified under DIN EN ISO 9001 and 13485.



Test mark for
quality management system



GREINER

GREINER GmbH
Wettestraße 1 · 74385 Pleidelsheim · Germany
T +49 7144 81120 · F +49 7144 8112 99
info@greiner-gmbh.de · www.greiner-gmbh.de